



Makers Making Change

A Neil Squire Program



LipSync Gaming SETUP GUIDE

Table of Contents

Table of Contents	2
About the LipSync Gaming	3
What Do You Get?	3
Types of LipSync	3
Device Compatibility	4
Powering the Controller	4
LipSync Gaming Firmware Installation	5
Setting Up the LipSync Gaming for PC Gaming	7
Functions (PC Games)	7
Setting Up the LipSync Gaming with the Xbox Adaptive Controller	8
Functions (Xbox Adaptive Controller)	10
Customize Profile	10
Cleaning the LipSync	13
References	13

PROGRAM OF:

FUNDING PROVIDED BY:



Foundation




About the LipSync Gaming

The LipSync Gaming emulates a gaming controller for PC gaming and the XBOX console. Users can move the mouthpiece up, down, left, and right; and use the sip-and-puff for other controller functions. This adaptive gaming device is suitable for users with limited or no hand and finger dexterity. It can be paired with auxiliary controllers, like the XBOX Adaptive Controller, to suit a gamer's needs.

It is an open source assistive technology design by Makers Making Change, a program of the Neil Squire Society that connects volunteer makers with people with disabilities to create more affordable assistive technology. Join the Makers Making Change to help people in your community: <http://www.makersmakingchange.com>


What Do You Get?

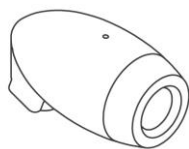

LipSync



Filter x 3

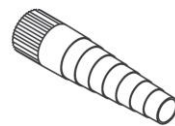

Mouthpiece x 3

Types of LipSync

	Description	Connection Type
LipSync	Works like a mouse	USB
LipSync Wireless	Works like a mouse	Wireless
LipSync Macro	Works like a keyboard	USB & wireless
LipSync Gaming	Works like a gaming joystick	USB

Device Compatibility



Windows
MacOS



Xbox
(only with the Xbox
Adaptive Controller)



Xbox Adaptive
Controller

Powering the Controller


The Adaptive Controller has an internal rechargeable battery that allows users to use it wirelessly. Using the LipSync (and any external accessibility switch) will draw additional power.

To prevent fast power consumption, we recommend the Adaptive Controller be plugged into to an electrical outlet using an AC Adapter (Option 1) or a power bank of at least 5000mAH (Option 2) using a USB to USB-C cable while gaming.

Option 1




Option 2

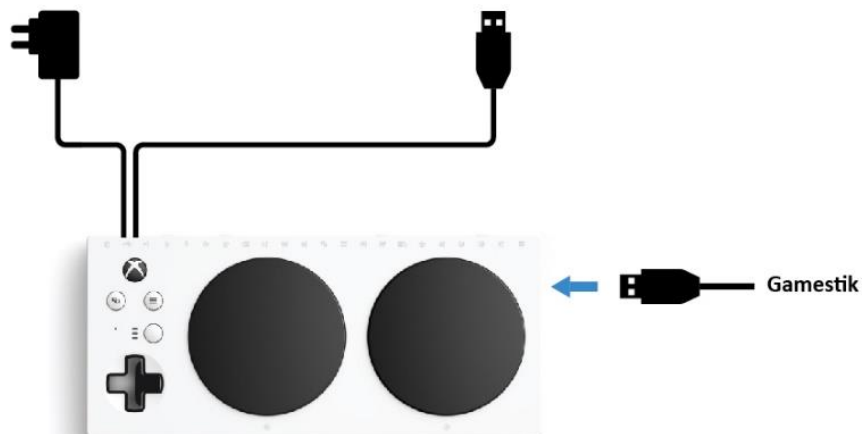

LipSync Gaming Firmware Installation

1. Download and install the latest version of Arduino IDE from the following link:
<https://www.arduino.cc/en/main/software>
2. Download the Gaming firmware from the following link:
https://github.com/makersmakingchange/LipSync-Gaming/raw/master/Software/LipSync_Gaming_Firmware/LipSync_Gaming_Firmware.ino
3. Visit the *"ArduinoJoystickLibrary"* GitHub repository using following link: [1]
<https://github.com/MHeironimus/ArduinoJoystickLibrary>
4. Download the *"ArduinoJoystickLibrary"* by clicking on *"Clone or Download"* button and then selecting *"Download ZIP"* option.
5. Extract the "ArduinoJoystickLibrary-master.zip" file.
6. Copy the Joystick folder to the Arduino libraries folder. An example of path to Arduino libraries folder in Windows: *"C:\Program Files (x86)\Arduino\libraries"*.
7. Once the folder is copied, the Joystick library should appear in the Arduino IDE list of libraries.
8. Plug the LipSync USB cable into the computer port.
9. Open Arduino IDE and you should be able to see Joystick examples under *"File>Examples>Joystick"*.
10. Now open the *"LipSync_Gaming_Firmware.ino"* file using Arduino IDE.
11. Select the *"Arduino/Genuino Micro"* as the Board and then select the available port number under Tools option.
12. Click on upload button to upload the code.

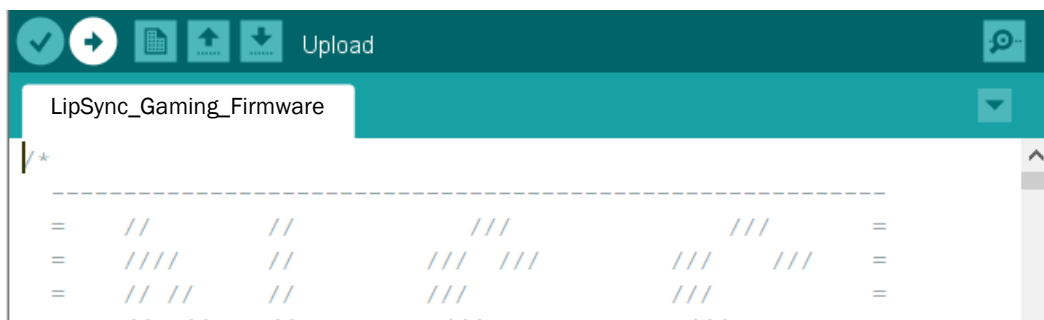


13. With the LipSync Gaming plugged into a power source, perform the joystick calibration.



Setting Up the LipSync Gaming for PC Gaming

1. Plug the LipSync USB to any computer or laptop USB Port.
2. Start up the PC game. The LipSync will not behave like a computer mouse outside of the gaming window.
3. Adjust the sensitivity of the sip-and-puff features by pressing the triangular buttons on the back of the LipSync. Increase the sensitivity by pressing the ▲ triangular button and decrease it by pressing the ▼ triangular button.

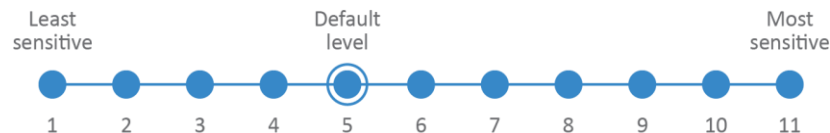

4. If any of the functions starts to drift while playing a game, release the mouthpiece and hold both ▲ ▼ buttons at the same time for 1 second to perform the joystick initialization.

Functions (PC Games):

LipSync Action	Time (Sec)	Mode 1	Mode 2
		(Hold ▼ for 2 sec) ¹	(Hold ▲ for 2 sec) ²
Mouthpiece (Right)	N/A	Move right	Move right
Mouthpiece (Down)	N/A	Move down	Move down
Mouthpiece (Left)	N/A	Move left	Move left
Mouthpiece (Up)	N/A	Move up	Move up
Puff	1	Button Number 1	Button Number 1
Sip	1	Button Number 2	Button Number 2
Puff	2	Button Number 3	N/A
Sip	2	Button Number 4 ³	N/A
Puff	3	Button Number 5	N/A
Sip	3	Button Number 6	N/A


¹ LED will blink red and green 2 times


² LED will blink red and green 1 time

³ The press button 1 or 2-second long sip is reserved to be used for Xbox Adaptive Controller.




Setting Up the LipSync Gaming with the Xbox Adaptive Controller

1. Turn on your XBOX Console.
2. Turn on your Adaptive Controller by pressing the XBOX button .
3. Plug the LipSync USB Cable to the right or left USB port on either side of the Adaptive Controller. **Each side has different gaming functions** so double check with the user.

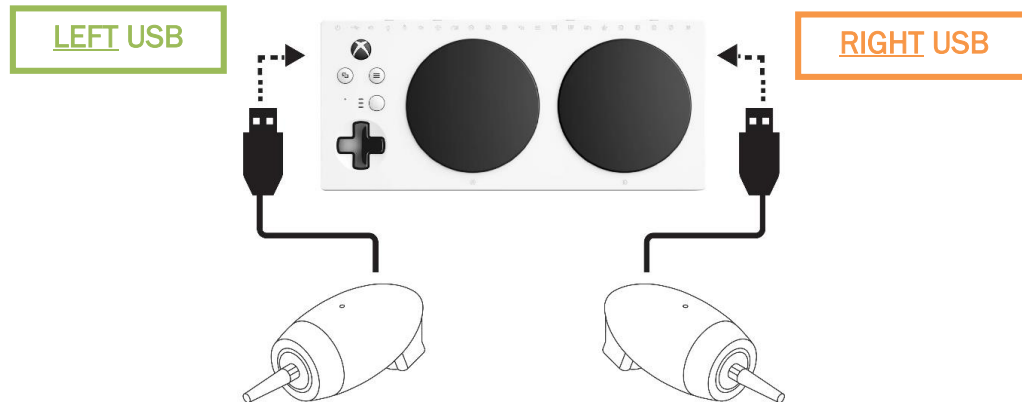

4. Adjust the sensitivity of the sip-and-puff features by pressing the triangular buttons on the back of the LipSync. Increase the sensitivity by pressing the ▲ triangular button and decrease it by pressing the ▼ triangular button.

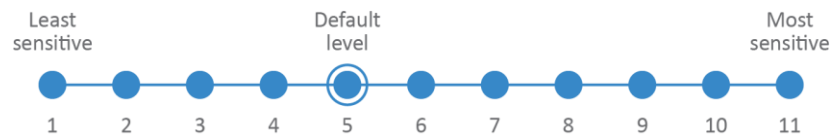

5. If any of the functions starts to drift while playing a game, release the mouthpiece and hold both ▲ ▼ buttons at the same time for 1 second to perform the joystick initialization.

Functions (Xbox Adaptive Controller):

The following Mode 2 functions are the default options for the specified actions which can be modified using profile settings.

Action	Time (Sec)	Mode 1 (Left USB Port) (Hold ▼ for 2 sec) ¹	Mode 2 (Left USB Port) (Hold ▲ for 2 sec) ²	Mode 1 (Right USB Port) (Hold ▼ for 2 sec) ¹	Mode 2 (Right USB Port) (Hold ▲ for 2 sec) ²
Right	N/A	Move right	Move right	Move right	Move right
Down	N/A	Move down	Move down	Move down	Move down
Left	N/A	Move left	Move left	Move left	Move left
Up	N/A	Move up	Move up	Move up	Move up
Puff	1	X1	X1	View	View
Sip	1	X2	X2	Menu	Menu
Puff	2	LS	N/A	RS	N/A
Sip	2	LB / Shift	N/A	RB	N/A
Puff	3	A	N/A	X	N/A
Sip	3	B	N/A	Y	N/A

¹ LED will blink red and green 2 times

² LED will blink red and green 1 time

Customize Profile

The Xbox Adaptive Controller comes with the option to configure and map buttons using 3 different mapping profiles. Switch between profiles by pressing on the Profile Button or pressing the **Shift** button after it's assigned in settings. [2]



Please follow the instructions below to map the buttons through the **Profile** settings:

1. Press the **Xbox** button to open the guide.
2. Select Home > Settings > Kinect and Devices > Devices & Accessories.
3. Press the **A** button (3 second puff) > Select **Xbox Adaptive Controller**.
4. Select **Configure** option > press the **A** button (1 second puff).
5. Select **New profile** > press the **A** button (3 second puff).
6. Enter the name for the new profile (Example: "**LipSync1**").
7. Under the **Mapping** menu, select the Xbox controller button that you would like to map > press the **A** button (3 second puff).
8. Assign the **Primary** and **Shift** button actions for selected button > press the **B** button (3 second sip).
9. The following five examples can be used to map the main buttons:
 - **X1** button or short puff
 - A. Under the **Mapping** menu, select the **X1** button > press the **A** button (3 second puff).
 - B. Assign the **Primary** and the **Shift** button actions for selected **X1** button.
 - C. The **Primary** button action is assigned to the **A** button and the **Shift** button action is assigned to the **X** button for the **X1** button in the **LipSync1** profile example.
 - D. Press the **B** button (3 second sip).

- **X2** button (1 second sip)
 - A. Under the **Mapping** menu, select the **X2** button > press the **A** button (3 second puff).
 - B. Assign the **Primary** and the **Shift** button actions for selected **X2** button.
 - C. The **Primary** button action is assigned to the **B** button and the **Shift** button action is assigned to the **Y** button for the **X2** button in the **LipSync1** profile example.
 - D. Press the **B** button (3 second sip).
- **LS** button or 2 seconds long puff
 - A. Under the **Mapping** menu, select the **LS** button > press the **A** button (1 second puff).
 - B. Assign the **Primary** and the **Shift** button actions for selected **LS** button.
 - C. The **Primary** button action is assigned to the **LS** button and the **Shift** button action is assigned to the **RS** button for the **LS** button in the **LipSync1** profile example.
 - D. Press the **B** button (3 second sip).
- **A** button or 3 seconds very long puff
 - A. Under the **Mapping** menu, select the **A** button > press the **A** button (1 second puff).
 - B. Assign the **Primary** and the **Shift** button actions for selected **A** button.
 - C. The **Primary** button action is assigned to the **View** button and the **Shift** button action is assigned to the **View** button for the **A** button in the **LipSync1** profile example.
 - D. Press the **B** button (3 second sip).
- **B** button or 3 seconds very long sip
 - A. Under the **Mapping** menu, select the **B** button > press the **A** button (1 second puff).
 - B. Assign the **Primary** and the **Shift** button actions for selected **B** button.
 - C. The **Primary** button action is assigned to the **Menu** button and the **Shift** button action is assigned to the **Menu** button for the **B** button in the **LipSync1** profile example.
 - D. Press the **B** button (3 second sip).

10. Under the **Mapping** menu, select the button that you would like to map to the **Shift** button action > press the **A** button (1 second puff).

11. Select the **Use as Shift** button check box > press the **A** button (1 second puff) > press the **B** button (1 second sip).
12. Press the **B** button (1 second sip) to enter the **Profiles** page again.
13. Select the new created **LipSync1** profile > press the **A** button (1 second puff).
14. Save the new profile by selecting the **No Slot** dropdown menu > press **the A** button (1 second puff).
15. Select one of the three available slots > press the **A** button (1 second puff).
16. The **LipSync1** profile is now ready to use. Press the Profile button on the Xbox Adaptive Controller to switch between profiles and select **LipSync1** profile for these functions:

Action	Time (Seconds)	Mode 1 Primary	Mode 1 Shift
Mouthpiece (Right)	N/A	Move right	Move right
Mouthpiece (Down)	N/A	Move down	Move down
Mouthpiece (Left)	N/A	Move left	Move left
Mouthpiece (Up)	N/A	Move up	Move up
Puff	1	A	X
Sip	1	B	Y
Puff	2	LS	RS
Sip	2	LB (Shift Button)	LB (Shift Button)
Puff	3	View	View
Sip	3	Menu	Menu

Cleaning the LipSync

The mouthpiece and filter are medical grade products used for capturing saliva. Each user should have their own set of mouthpiece and filter for hygienic reasons. It has a life cycle of 3-6 months depending on how often it's used and how much saliva is trapped by the filter.

Unused pieces should be kept in a clean sealed bag away from the sun and used ones should be disposed of after each cycle.

The GamesStik casing is 3D printed layer by layer in plastic. Due to this, tiny gaps are formed between each layer, making the casing porous and not resistant to water.

Do:

- Wipe with a non-abrasive cloth lightly dampened with water
- Benzalkonium chloride antiseptic wipes can be used on mouthpiece
- Mouthpiece can be rinsed in cold water and air dried in its 3-6 month life cycle

Don't:

- Avoid using abrasive cleaning materials
- Alcohol will destroy the plastic material
- Do not use any cleaning agents or soap
- Do not wash the filter

References

1. [ArduinoJoystickLibrary Github repository](#)
2. [Customize the Xbox Adaptive Controller in the Xbox Accessories app](#)
3. [How do I use Shift with the Xbox Adaptive Controller?](#)